


Fast canine-retraction during levelling phase



The Hybrid Retractor achieves true bodily movement by applying a rotational and up-righting force. The superelastic spring in the Hybrid Retractor applies a constant physiologically correct force on the canine. The Hybrid Retractor is suitable for treatment in the upper and lower jaw and is used in combination with a vertical slot (.018 x .018) bracket for anchorage.

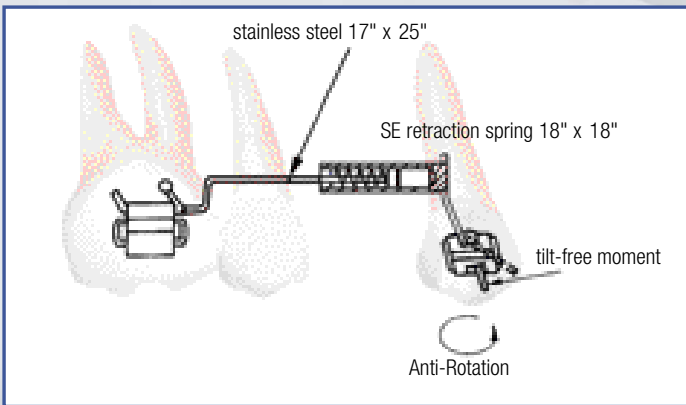

Hybrid Retractor: Fast retraction of canines using physiologically correct force with a superelastic spring avoiding unwanted tipping and rotation.

Retraction force: approx. 100-150 g
Spring load: 5 mm

Hybrid retractor*



Hybrid retractor € 0297

Order Nr.	content	description
307-1020	2	Hybrid retractor



€ 0297

Maxillary-kits € 0297

	Order Nr.	tooth	content
.018"	307-1030	<u>3</u> / <u>3</u>	2 Brackets for bonding horizontal slot .018", vertical slot .018" x .018" und 2 Hybrid retractors
.022"	307-1035	<u>3</u> / <u>3</u>	2 Brackets for bonding horizontal slot .022", vertical slot .018" x .018" und 2 Hybrid retractors

Mandibular-kits € 0297

	Order Nr.	tooth	content
.018"	307-1040	<u>3</u> / <u>3</u>	2 Brackets for bonding horizontal slot .018", vertical slot .018" x .018" und 2 Hybrid retractors
.022"	307-1045	<u>3</u> / <u>3</u>	2 Brackets for bonding horizontal slot .022", vertical slot .018" x .018" und 2 Hybrid retractors



Bernhard Förster GmbH · Westliche Karl-Friedrich-Straße 151 · 75172 Pforzheim · Germany
 Phone + 49 7231 459-0 · Fax + 49 7231 459-102 · info@forestadent.com · www.forestadent.com
 Forestadent USA · 2301 Weldon Parkway · St. Louis, MO 63146 · TOLL FREE 1-800-721-4940
 Fax 314-878-7604 · info@forestadentusa.com · Internet: http://www.forestadentusa.com
 Forestadent Ltd. · 21 Carters Lane/Kiln Farm · Milton Keynes MK11 3HL · Great Britain
 Phone + 44 19 08 56 89 22 · Fax + 44 1 908 560 611 · info@forestadent.co.uk

*acc. to Prof. Dr. F. G. Sander
 Infosheet No. 171/9.2005