




*Ceramic doesn't
have to splinter!*

*Do you have problems when removing
ceramic brackets? If you do read on ...*




Quicklear®-Brackets

Patents no: EP 0782414 - US 5820371




Quicklear®-brackets are easily and quickly removed without fracturing

Have you ever experienced the problem of ceramic brackets Breaking or fracturing when removing them and the time-consuming task of trimming away the residue? This is now a thing of the past. Try the new Quicklear® ceramic brackets and discover easy, splinter-free removal, which is 99.9% guaranteed by FORESTADENT.

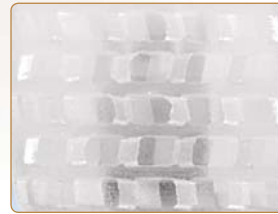

Quicklear® ceramic brackets have a newly developed, inversive bracket base with diametrical rows of holes. The patent-protected base ensures not only excellent bond strength but also easy removal of the bracket upon completion of treatment.

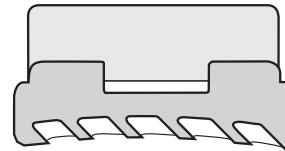

Enamel chipping is prevented, as a large amount of adhesive remains on the tooth after removal of the bracket.

Quicklear®-brackets are removed this easily...




Hold the new specially designed "Pauls Tool" parallel to the occlusal plane ...



... and place the vertical bars mesiodistally to the bracket body.



Push the end of the "Pauls Tool" towards the distal and remove the bracket. The entire bracket is removed by the new debonding instrument without fracturing, ...

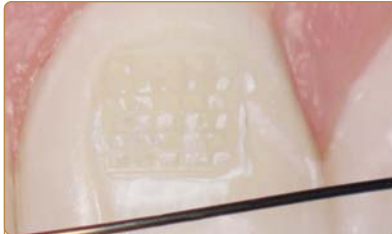

... while the adhesive remains on the tooth. The adhesive can then be easily removed using an adhesive bur or cutter.



Brackets can be easily removed with this technique even when an archwire is inserted.



Alternatively, the bracket removing plier by FORESTADENT can be used for debonding.

Features and benefits of *QuicKlear*[®]-brackets



The QuicKlear[®] bracket was developed as a self-ligating, translucent ceramic bracket and is an excellent option in addition to the QuicKlear[®]2.0[®] bracket system. The outstanding aesthetics and stability of the QuicKlear[®] bracket ensure a natural smile and short chairside visits.

The following features guarantee doctor and patient satisfaction:

- outstanding translucency
- highly efficient force transfer
- purely mechanical retention provided by the patented bracket bases
- easy splinter-free debonding of the brackets upon completion of treatment
- four rounded contact points in the slot prevent binding and notching
- rounded bracket and slot edges ensure intraoral comfort and minimal friction
- interactive locking clip can be easily opened and closed

For more information about QuicKlear[®] please contact your Forestadent USA sales representative or call 1-800-721-4940 for more details.



Are you interested in solutions without concerns?

Try **QuicKlear**[®]-bracket starter-Kits

Introducing the new starter kits.



QuicKlear/Quick 2.0 starter kits – 3 cases

options:

QuicKlear Roth .018" 3-5 with hooks/Lower Quick 3 with hooks or
QuicKlear Roth .018" 3-5 with hooks/Lower Quick 3-5 with hooks or
QuicKlear Roth .022" 3-5 with hooks/Lower Quick 3 with hooks or
QuicKlear Roth .022" 3-5 with hooks/Lower Quick 3-5 with hooks

All kits contain:

- 3-case sets of QuicKlear/ Quick brackets
- BioStarter upper archwire .012" 5 archwires
- BioStarter lower archwire .012" 5 archwires
- QuicKlear model
- QuicKlear & Quick special probe
- QuicKlear bracket debonder "Pauls-Tool"
- Archwire director
- Clamping tube open
- Translucent QuicKlear model

Price per kit:
\$ 1,280.00

QuicKlear/Quick 2.0 starter kits – 6 cases

options:

QuicKlear Roth .018" 3-5 with hooks/Lower Quick 3 with hooks or
QuicKlear Roth .018" 3-5 with hooks/Lower Quick 3-5 with hooks or
QuicKlear Roth .022" 3-5 with hooks/Lower Quick 3 with hooks or
QuicKlear Roth .022" 3-5 with hooks/Lower Quick 3-5 with hooks

All kits contain:

- 6-case sets of QuicKlear/ Quick brackets
- BioStarter upper archwire .012" 6 archwires
- BioStarter lower archwire .012" 6 archwires
- QuicKlear model
- QuicKlear & Quick special probe
- QuicKlear bracket debonder "Pauls-Tool"
- Archwire director
- Clamping tube open
- Translucent QuicKlear model

Price per kit:
\$ 1,520.00

If our starter kits are too expensive, try our smaller kit option!

QuicKlear[®]-bracket introductory offer

QuicKlear bracket 3-case kit
upper 5 to 5 incl. "Pauls Tool"

3-case-kit:
\$ 550.00



TOLL FREE Order: 1-800-721-4940



GERMAN PRECISION IN ORTHODONTICS

2315 Weldon Parkway · St. Louis, MO 63146 · **TOLL FREE 1-800-721-4940**

Fax 314-878-7604 · info@forestadentusa.com · www.forestadentusa.com